




Please read the following QR code for the Japanese and Chinese manuals for the "eyenee" product.



Click here for product information!

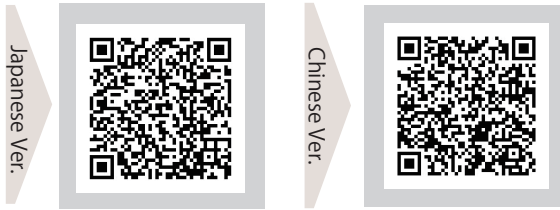

CRECHES Co., Ltd.  
1-6-7 Azuchimachi, Chuo-ku, Osaka 541-0052

## Instruction Manual



Before using this product, please read this instruction manual carefully and use

## Getting Started

Thank you for purchasing the eyenee.

To use this product correctly, be sure to read this instruction manual before use.

Also, after reading the manual, a warranty card is attached at the end of the manual.

Please keep it in a safe place so that it will not be lost or soiled.

### Table of contents

1.Safety precautions	P1~3
2.Explanation of this device	P4
3.Mode function of this device	P5
4.How to charge this device	P5~6
5.Procedure for using this device	P6~7
6.Treatment method	P7
7.Handling of adhesive pads	P8
8.Product specifications	P9
eyenee warranty	P10

### 1.Safety Precautions


- Please read this "Safety Precautions" carefully before use and follow the instructions for proper use.
- The precautions given here are for the safe and correct use of the product and for prevention of harm and damage to the user and others.

#### The notes are classified as follows.

These are all important safety items, so be sure to follow them.

Example of picture display

- ⚠ Danger: Contents that are expected to cause imminent risk of death or serious injury to a person.
- ⚠ Warning: Content that may lead to death or serious injury.
- ⚠ Caution: Contents that assume the risk of human injury and the occurrence of physical damage.
- 🚫 Prohibition: This is a notification of a prohibited act.


⚠ Danger	<p>🚫 <b>Do not use with the following medical electronic devices. Be sure to remove precious metals and watches.</b></p> <p>① Implantable medical electronic devices such as pacemakers ② Medical electronic devices for life support such as heart-lung machines ③ Wearable medical electronic devices such as electrocardiographs</p> 
	<p>🚫 <b>The following people should avoid using this device.</b></p> <p>① People with heart disease ② People who are pregnant or may be pregnant ③ People who are drinking alcohol</p>
	<p>🚫 <b>For usage time and frequency, please follow the directions in the instruction manual.</b></p> <p>⚠ Excessive use may cause physical injury.</p>
⚠ Warning	<p>🚫 <b>The following people, please be sure to consult with your doctor before use.</b></p> <p>① Hospital patients ② People with atomic-bomb sickness ③ People with malignant tumor ④ People with acute illness ⑤ People with abnormal blood pressure ⑥ 10 years old or younger ⑦ People with infectious disease ⑧ People with tuberculosis ⑨ People with hemophilia ⑩ People with legal infectious disease ⑪ People with melasma facialis ⑫ People with allergic diseases ⑬ People taking medication ⑭ People with extreme fatigue ⑮ People with febrile illness having temperature of 38°C or higher ⑯ People who have abnormal skin ⑰ People with capillary dilation due to prolonged use of steroid hormone drugs or impaired liver function ⑱ People with eye diseases ⑲ People who have undergone eye surgery or receiving eye treatment ⑳ Excluding LASIK ㉑ People who have capillary dilatation ㉒ People with cataract, glaucoma, retinal detachment ㉓ If you feel any abnormality during use, please stop using immediately.</p>
	<p>🚫 <b>Please do not attach the adhesive pad close to the heart. The heartbeat may get disrupted.</b></p>
	<p>🚫 <b>Avoid using on the following body parts.</b></p> <p>① Mucous membrane part ② Anus and genital region ③ Abdominal part during menstruation ④ Plastic surgery part ⑤ Eyeball ⑥ Joints such as knees and elbows ⑦ Wound or site just after hair loss ⑧ Near the heart ⑨ Near thyroid gland, head, or inside the mouth ⑩ Parts with pain/sensory disturbances ⑪ Parts with acne and pimples causing inflammation ⑫ Parts with skin inflammation due to cosmetics ⑬ Parts with atopic/eczema symptoms ⑭ Itching and burning, parts with pathological spots due to physical irritation, etc. ⑮ Parts with metal, plastic, silicon, etc. embedded in the body</p>
⚠ Warning	<p>🚫 <b>Do not repair, disassemble, or modify the main unit and accessories without permission.</b></p> <p>⚠ You may not be able to receive maintenance or repair even during the warranty period.</p>
	<p>🚫 <b>Do not use for any purpose other than its original purpose.</b></p> <p>⚠ It may cause an accident or malfunction.</p>
	<p>🚫 <b>If each pad or cord is broken or if any abnormality occurs on the main unit, immediately stop using the product and be sure to request repair.</b></p> <p>⚠ It may cause fire, burns or accidents.</p>
	<p>🚫 <b>Do not use with other treatment devices at the same time, or use with an application agent (including spray type).</b></p> <p>⚠ You may feel sick or get sick.</p>
	<p>🚫 <b>Do not use this product with other electrical products, or use or store it in a place where strong radio waves or magnetism is generated.</b></p> <p>⚠ It may cause malfunction or failure due to interference.</p>
	<p>🚫 <b>Do not remove the adhesive pad during treatment. Please make sure to pause the treatment once when you take it off and on.</b></p> <p>⚠ You may receive a strong shock.</p>
<p>🚫 <b>Do not use near highly volatile substances (benzine, thinner, nail polish remover).</b></p> <p>⚠ There is a risk of explosion and fire.</p>	

- Prohibition** ● When using the device in the bath, please do not soak the main unit in hot or cold water.  
▼ Leaving the device in a humid place such as bathroom, etc. for several days may cause a malfunction.
- Prohibition** ● Do not let any metal objects get in contact with the adhesive pad and EMS electrode.  
▼ You may receive a strong shock.
- Prohibition** ● Do not use the product while driving a car.  
▼ May cause strong irritation, accidents or malfunctions.
- Prohibition** ● Please do not use while sleeping.  
▼ There is a risk of wire disconnection accidents and troubles.
- Prohibition** ● If you notice any gas leak, immediately stop using the product, ventilate the room, turn off the power, and immediately disconnect the AC adapter from the outlet.  
▼ The gas may catch fire and cause an explosion or fire.
- Prohibition** ● If you hear thunder, stop using the product immediately.  
▼ Lightning strikes may cause an overcurrent and abnormal output.
- Prohibition** ● Do not use or leave the product in places with high temperatures (such as places exposed to direct sunlight or near heating devices such as stoves), places with a lot of humidity or dust, or places exposed to steam (such as a cooktop, bathroom, or humidifier).  
▼ It may cause fire or electric shock.
- Prohibition** ● Do not place the product on unstable places, drop, or place objects on the product. Also, do not push or hit the operation section more than necessary.  
▼ It may cause malfunction or injury.
- Prohibition** ● Please keep out of reach of children.  
▼ There is a risk of injury.


- Prohibition** ● When removing the charging USB cord or adhesive pad cord, do not hold it with wet hands or pull it forcibly.  
▼ You may get an electric shock or a short circuit may cause a fire.
- Prohibition** ● Do not use if the charging USB cord is damaged or loose when inserted into electrical outlet. Also, do not use with multi-socket adapters.  
▼ You may get an electric shock or a short circuit may cause a fire.
- Prohibition** ● Do not damage, break, modify, forcibly twist, bend, stretch, or tie up the cords in a bundle.  
▼ There is a risk of electric shock and fire.
- Prohibition** ● Although this unit is waterproof to the extent that it may tackle wet conditions, rain, etc., please note that it cannot be used underwater or in the sea.  
▼ You may get an electric shock or a short circuit may cause a fire.
- Prohibition** ● Disconnect the USB charger cord from the main unit unless during charging.  
▼ Injury, burns or insulation deterioration may cause electric shock, short circuit, or fire.

## 2.Explanation of this device

### ●Standard accessories

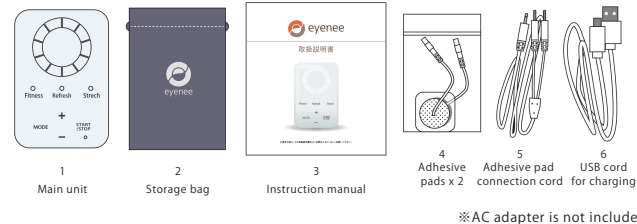

### ●Name of each part

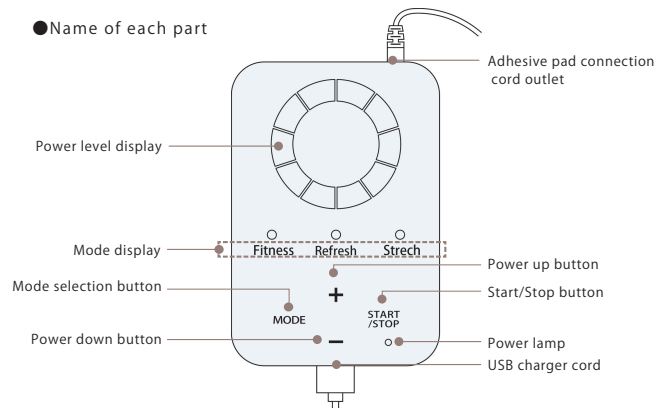

### 3. Mode function of this unit

Mode	Contents	Time
Fitness	EMS Fitness system helps to exercise and train the muscles around the eyes.	30 minutes
Refresh	EMS Refresh system refreshes the muscles around the eyes.	30 minutes
Stretch	EMS Stretching system comfortably stimulates the muscles around the eyes.	30 minutes

The timer for each mode is pre-configured and cannot be set arbitrarily.

### 4. How to charge the main unit


This device cannot be used while being charged.

After charging is complete, remove the charging cable before use.

If you use this product for the first time after purchase or if you have not used it for a long time, charge it for about 2.5 hours.

#### (Guidelines for charging)

- When the remaining charge is 5% or less, it takes about 2.5 hours to fully charge the battery.
- You can use the device for 10 hours or more with one full charge.

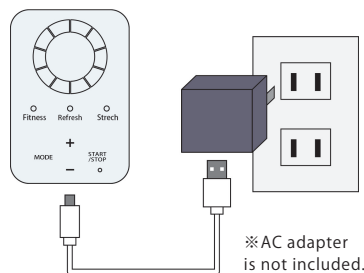

- 1 Insert the USB charger cord into the charging connector on the bottom of the main unit.
- 2 Insert the USB charger cord into your AC adapter to charge.
- 3 While charging, all the power level indicators light up and the power lamp lights up.
- 4 When charging is complete, the power level indicator turns off.


#### Low power display

- When the remaining power is 20% or more, the power lamp is on.
- When the remaining power is 20% or less, the power lamp blinks every 1.5 seconds.
- When the remaining power is 5% or less, the power lamp blinks at 0.5 second
- Room temperature suitable for charging is 15 to 30°C. Charging at low temperatures can shorten the battery life, so charge at an appropriate temperature. If it is below 0°C, the charging speed may decrease.
- To protect the built-in battery, charge the battery once a month even if it hasn't been used for a long time.
- If the operating time is short even when the battery is properly charged, it means that the battery is dead.

#### After charging

Unplug the AC adapter from the outlet and disconnect the USB charger cord from the charging connector.

### 5. Procedure for using this device

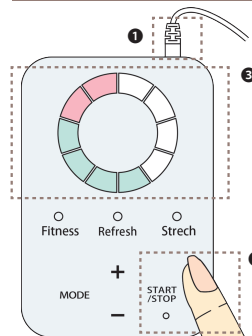


#### Connecting the adhesive pad connection cord

- 1 Insert the adhesive pad connection cord into the jack outlet on the top of ① and connect it.

#### How to turn on/off the power

- 2 If you press and hold START / STOP on ②, a beep will sound twice, the circular lamp in ③ will go around once, and the power will turn on.
- 3 Press and hold START / STOP again to turn off the power.
- 4 If no operation is performed for 15 minutes, the power will turn off automatically.

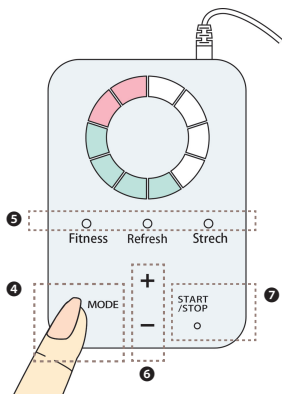


### Mode selection method

- 5 When you press MODE on ④, the mode indicator lights up in the order of [Fitness] → [Refresh] → [Stretch].
- 6 Select one of the three modes and if you press START / STOP on ⑦ the selected mode will change the lamp indicator from steady lit to blinking. ※Once the mode starts, you cannot switch to other modes.
- 7 Press + and - of ⑥ to adjust the power. ※You can also adjust the power by pressing + or - during use.
- 8 After 30 minutes, a beep sounds and the operation ends.
- 9 If you want to stop during treatment, press START / STOP in ⑦ will change the lamp indicator from blinking to steady lit and then you can stop.


### 6.Treatment method

Firmly insert the "adhesive pad terminal" and the "adhesive pad connection cord" and attach the adhesive pad to the treatment area. Please perform treatment at a comfortable level where you can feel the contraction of muscles caused by electricity.

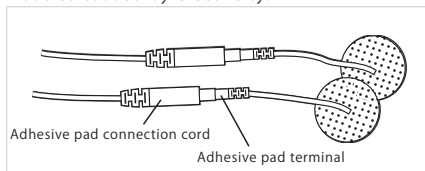


### 7.Handling of adhesive pads

#### Precautions for pasting adhesive pads

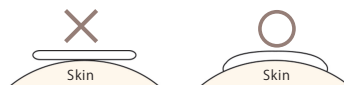

Make sure the adhesive pad is properly glued to your skin.

#### How to store adhesive pads

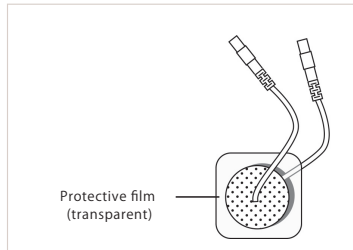

- After use, be sure to put a protective film (transparent) on the adhesive surface to preserve it. If the product is stored without the protective film attached, the adhesive surface is sticky, so dirt and dust are likely to adhere, and the adhesive strength may be weakened, and the product life may be shortened.
- Do not scratch the adhesive surface with your fingernail, brush, etc.
- Do not use hot water, detergent, chemicals (thinner, benzene, etc.).

- Danger**
- Please do not attach the adhesive pad close to the heart. The heartbeat may get disrupted.
  - If you use the product with the adhesive surface damaged or peeled off, it may cause burn injuries. Be sure to replace it with a new one before use.

- Caution**
- Remove precious metals such as necklaces, piercings, rings, and watches before treatment.
  - Adjust the level according to your physical condition.
  - Be sure to use the special adhesive pad.
  - Make sure the adhesive pad is properly attached. You may feel a strong stimulus if there's an opening between the skin and the pad.
  - Set the output level low at the beginning, and gradually increase it as you get used to it. Do not suddenly do treatment at a high level.
  - Please perform treatment according to the condition of your skin. If you feel any abnormality on your skin, stop using it immediately.

## 8. Product specifications

### Main unit

Power supply	Rechargeable
Battery used	Lithium-ion battery 3.7V
Voltage	4.2V (MAX)
Powerconsumption	0.2W (MAX)
Unit size	H9.5×W6.5×D1.6cm
Weight	About 115g
Charging time	About 2.5 hours
Usage time	10 hours or more

## Warranty

This manual promise to carry out free repairs according to the warranty regulations.

### [Warranty Regulations]

1. If there is no description of customer name, store stamp, serial number, and date of purchase, the warranty will be invalid.
2. In the following cases, the product will be repaired for a fee even during the warranty period.
  - (A) Malfunction caused by incorrect use and improper repair or modification
  - (B) Malfunction and damage due to transportation, moving, drop, etc. after purchase
  - (C) When it is necessary to replace consumable parts
  - (D) If the warranty is lost, or the necessary information is unfilled
  - (E) When the content of the warranty is modified without permission
  - (F) Malfunction due to fire, earthquake, water damage, abnormal voltage, etc.
  - (G) When used for other than general household use (business use, etc.)
3. Please present this warranty to the retailer for repair. If you are unable to present the warranty due to loss, you will be charged for repair.
4. This warranty will not be reissued, so please keep it in a safe place.

Product name	eyenee (アイニー)
Serial number	
Purchase date	Year    Month    Date
Parts covered by warranty	1 year from the date of purchase•Main unit
Customer name	
Address	〒  TEL    (    )
Sales outlet	

Selling Agency    CRECHES Co., Ltd.  
1-6-7 Azuchimachi, Chuo-ku, Osaka 541-0052